

ClearHemograsper



Jaw design and function optimized for hemostasis

○ Diamond-shaped non-slip design

The tapered tip-type forceps design allows for easy grip without slipping.

○ 2 - Wires structure

By operating the jaw with a 2 wire and transmitting the gripping force to the end of the jaw, it helps for precise access to blood vessels without slipping.

○ Swing jaw

Easy access to lesions placed in anatomically difficult locations

○ Efficient heat transfer

Easy and accurate targeting to bleeding site


○ 1:1 rotation with the handle

Accurate bleeding site positioning and easy targeting

○ Excellent suction function & Irrigation port

Φ2.3 Excellent suction function with a thin sheath

Convenient when using scope without water jet function by inserting irrigation port

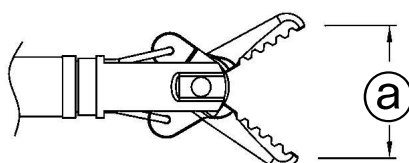

STERILE



Specification

| Model | Opening Width(@) | Feature | Working Length |
|-----------------|------------------|------------|----------------|
| FM-HG0001(1650) | 5mm | Rotational | 1,650mm |
| FM-HG0001(1820) | | | 1,820mm |
| FM-HG0001(2000) | | | 2,000mm |

Packaging : 5ea / box