



FineJet Cannula

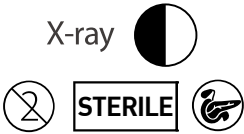

- 2 models available for the user's convenience
- Tube & Distal Tip design optimized for solution injection
- 4.4Fr Tapered Tip Design can be used for stenosis and small papilla, easy intubation
- Rounded Tip Forming
Minimize tissue damage
- Catheter tube with adequate stiffness
- Easy open/close handle design
Convenience when inserting a guide wire and injecting a contrast agent with a backflow prevention structure


Specification

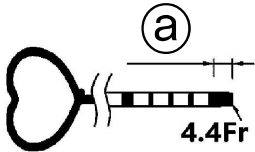

Model	Stylet	Tip Length(a)	Sheath Length	Tip O.D	Sheath O.D
FM-CN2540-R	Removable	5mm	2,200mm	Φ1.45mm	Φ1.7mm
FM-CN2540	Non-Removable				


Packaging : 5ea / box

Color Marker

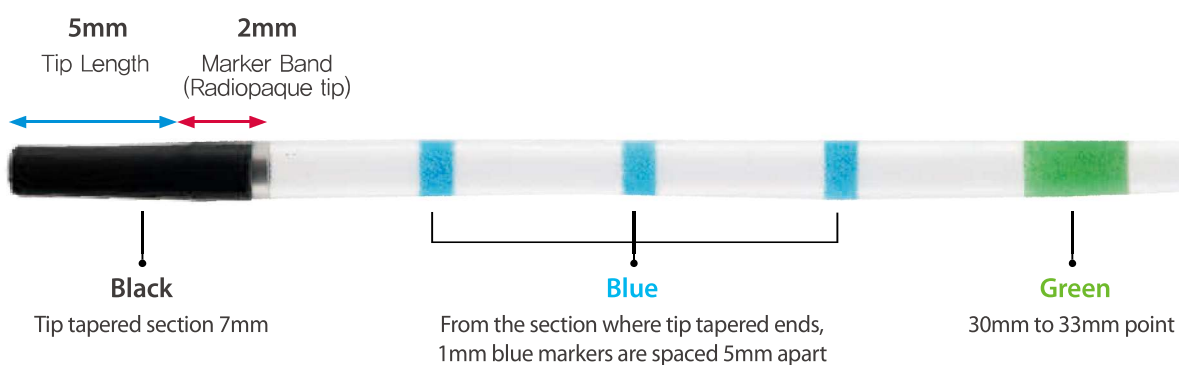

Open

Guidewire or Stylet can be freely inserted by pulling the slider

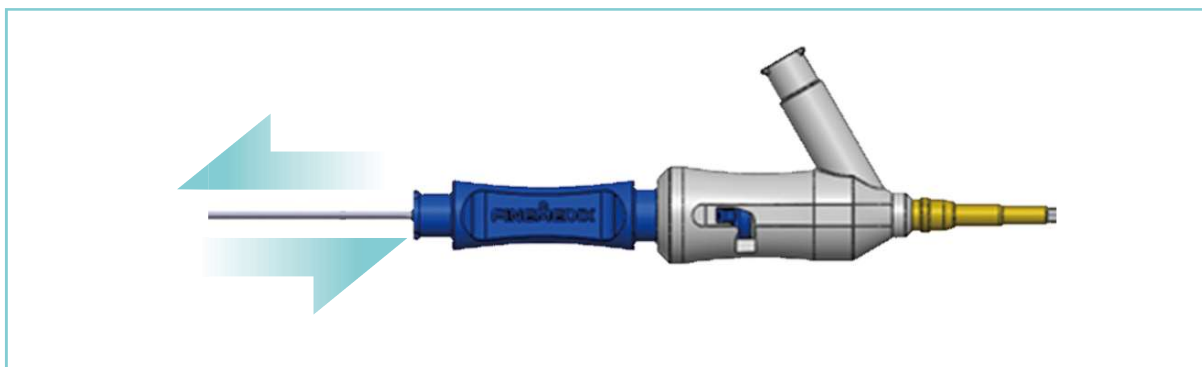

Close

By pushing the slider to lock it, can prevent backflow when the contrast medium is injected.

